

ELLAVI drape - detecting PPH

Attention: Place after delivery of baby and uterotonic administration, but before placenta delivery

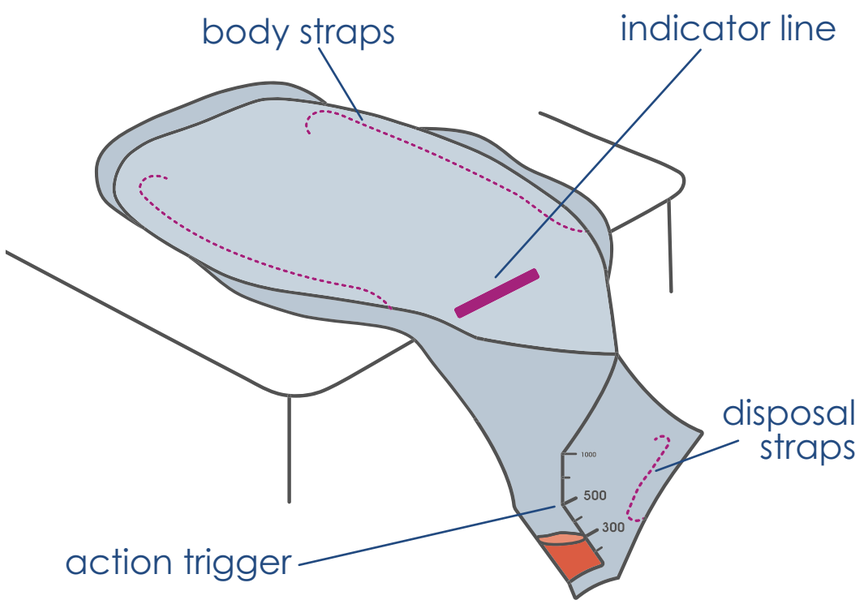


DISPENSER

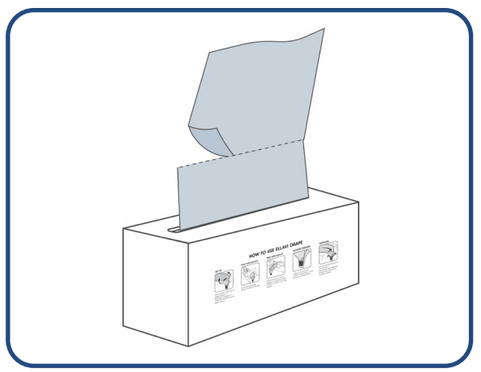


Remove drape from dispenser

SET UP

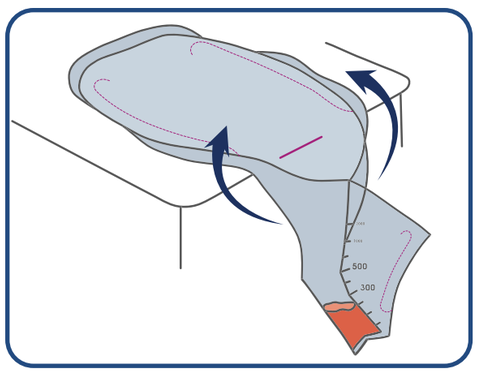


1) Place drape with line on edge of bed
2) Lift sides of drape to channel blood flow and open the pouch

FOOT END POUCH

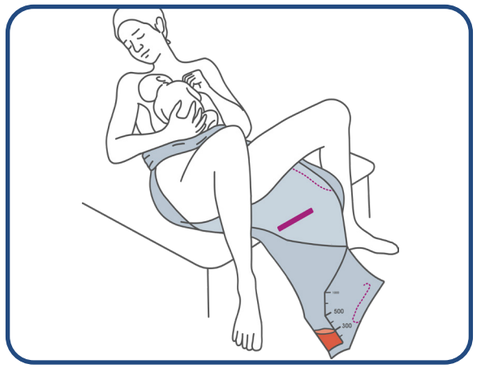


Tie body straps below the mother's breasts

BED SIDE POUCH

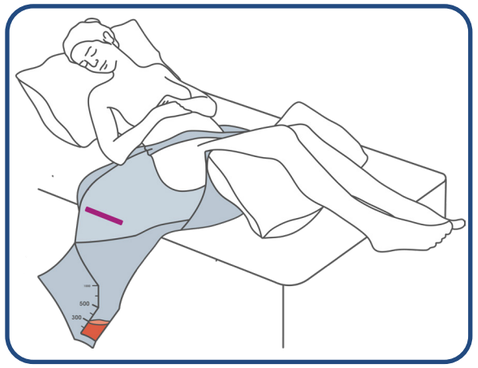

Tie body straps diagonally across the mother's abdomen


ACTION TRIGGER

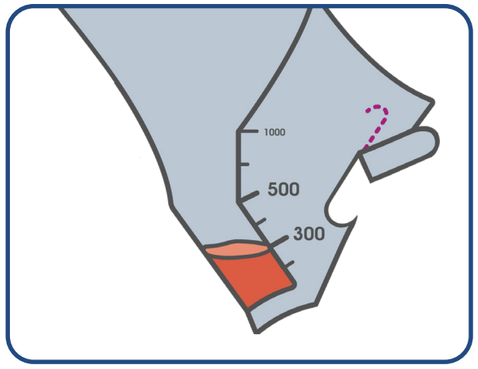

>500ml Blood loss

>300ml + one abnormal observation

DISPOSAL

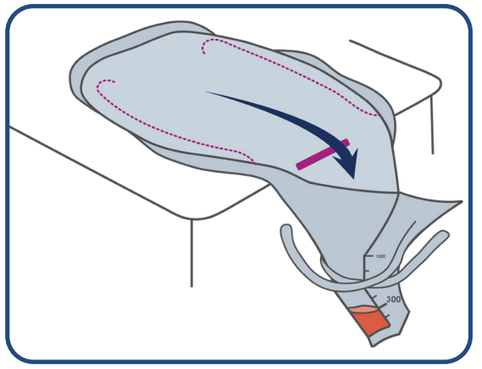

1) Scoop remaining blood into pouch
2) Tie pouch with disposal straps
3) Discard as medical waste

CAUTION: Federal law restricts this device to sale by or on the order of a physician.